

Dynabeads® Biotin Binder

Catalog no. 11047

Store at 2°C to 8°C

Publication No. MAN0007702

Rev. Date: January 2013 (Rev. 004)

Product Contents

Cat. no	Volume
11047	5 mL

Maximum product capacity

MNC*: $\sim 2 \times 10^9$ cells

Whole blood/buffy coat: ~ 200 mL

* Note: If using the product for negative isolation of multiple cell types simultaneously, the bead volume used is higher, thus giving a lower product capacity [see Table 1 and 2].

Dynabeads® Biotin Binder contains 4×10^8 beads/mL in phosphate buffered saline (PBS), pH 7.4, containing 0.1% bovine serum albumin (BSA) and 0.02% sodium azide as a preservative.

Caution: Sodium azide may react with lead and copper plumbing to form highly explosive metal azides.

Product Description

Dynabeads® Biotin Binder in combination with primary biotinylated antibodies are ideal for depletion or positive isolation of cells from different species (e.g. mouse, human), depending on the specificity of the primary antibody. Cells can be directly isolated from any sample such as whole blood, bone marrow, MNC suspensions or tissue digests. Other biotinylated molecules (e.g. peptides/proteins, lectins or nucleic acids) may also be used depending on the target.

The primary biotinylated antibodies are either added to the cell sample (indirect technique) or pre-coated onto the beads (direct technique) prior to cell isolation (fig. 1). The beads are then mixed with the cell sample in a tube. The beads bind to the target cells during a short incubation, allowing the bead-bound cells to be separated using a magnet.

Positive isolation

Discard the supernatant and use the bead-bound cells for downstream applications (e.g. isolation of proteins or nucleic acids [NA] or cell culture).

Note: If bead-free cells are required for downstream cellular applications or use in flow cytometry, use the equivalent CELLection™ Biotin Binder Kit (cleavage of a DNA-linker releases the cells from the beads). For isolation of phagocytic cells use Dynabeads® FlowComp™ Flexi (can be

used at the lower temperatures necessary to reduce phagocytic activity).

Depletion/negative isolation

Discard the bead-bound cells and use the remaining bead-free and untouched cells for applications. If different biotinylated monoclonal antibodies (mAb) are used to deplete several cell types simultaneously (negative isolation) to obtain untouched cells, use the indirect technique.

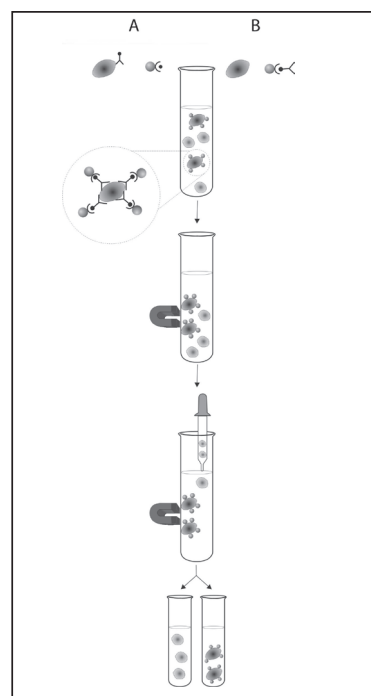


Figure 1: Cell isolation using (A) indirect or (B) direct technique.

Indirect versus Direct Technique

Use the *indirect technique* when depleting several cell types simultaneously with a cocktail of biotinylated mAbs (e.g. removing erythrocytes, platelets, and granulocytes from a MNC sample); very high depletion efficiency is required, the affinity of biotinylated antibodies is low, the cells have low expression of target antigens, or the direct technique gives unsatisfactory purity.

Use the *direct technique* when the primary antibody affinity is high, the cells have high expression of target antigens, or to make a larger stock preparation of primary coated magnetic beads (shelf life will be similar to that stated on the Dynabeads® vial).

Required Materials

- DynaMag™ Magnet (See www.lifetechnologies.com/magnets for recommendations).
- Mixer allowing tilting and rotation of tubes (e.g. HulaMixer® Sample Mixer).
- Isolation Buffer: Ca²⁺ and Mg²⁺ free PBS supplemented with 0.1% bovine serum albumin (BSA) and 2 mM EDTA, pH 7.4.
Note: BSA can be replaced by human serum albumin (HSA) or fetal calf serum (FCS). EDTA can be replaced by sodium citrate.
- Biotinylated antibodies.

General Guidelines

- Visit www.lifetechnologies.com/samplepreparation for recommended sample preparation procedures.
- Use a mixer that provides tilting and rotation of the tubes to ensure that Dynabeads® magnetic beads do not settle at the bottom of the tube.
- The choice of primary antibody is the most important factor for successful cell isolation. Note that some antibodies may show reduced antigen-binding efficiency when coated onto beads, even though the antibody shows good results in other immunological assays.
- To avoid non-specific binding of cells (e.g. monocytes, B cells), add aggregated IgG to block Fc-receptors prior to adding the primary antibodies.
- Wash cells prior to adding biotinylated antibodies or Dynabeads® magnetic beads to remove density gradient media (e.g. Ficoll) or soluble factors in serum (e.g. antibodies or cell surface antigens), which can interfere with the cell isolation protocol.
- This product should not be used with the MPC™-1 magnet (Cat. no. 12001D).
- Avoid air bubbles (foaming) during pipetting.
- Never use less than recommended volume of Dynabeads® magnetic beads.
- Carefully follow the recommended volumes and incubation times.
- Keep all buffers cold.

Protocols

Wash Dynabeads® Magnetic Beads

See Table 1 and 2 for volume recommendations.

- Resuspend the magnetic beads in the vial (vortex for >30 sec, or tilt and rotate for 5 min).
- Transfer the desired volume of magnetic beads to a tube.
- Add the same volume of Isolation Buffer, or at least 1 mL, and resuspend.
- Place the tube in a magnet for 1 min and discard the supernatant.
- Remove the tube from the magnet and resuspend the washed magnetic beads in the same volume of Isolation Buffer as the initial volume of magnetic beads (step 2).

Prepare Sample

- Cells can be directly isolated from any samples such as whole blood, bone marrow, MNC suspensions, or tissue digests. See "General guidelines" for sample preparation
- Prepare a MNC suspension according to "General Guidelines". Resuspend the cells at 1×10^7 cells/mL in Isolation Buffer.

Isolate Cells - Indirect Technique (Labeling Cells with Biotinylated Antibodies)

- See "General Guidelines" for recommendation of when to use the direct vs. indirect cell isolation technique.
- Use approximately 10 µg of primary biotinylated antibody per 10^7 target cells.

This protocol is based on 1×10^7 MNC or 1 mL whole blood, but is directly scalable from 1×10^7 to 4×10^8 cells or 1–40 mL whole blood. When working with lower volumes than 1×10^7 cells or 1 mL blood, use the same volumes as for 1×10^7 cells or 1 mL blood. When working with larger volumes, scale up all reagent and volumes accordingly, as shown in Table 1.

- Add ~ 10 µg primary antibody to 1 mL cell suspension and mix (titrate the antibody amount for your use).
- Incubate for 10 min at 2°C to 8°C.
- Wash the cells by adding 2 mL Isolation Buffer and centrifuge at $350 \times g$ for 8 min. Discard the supernatant.
- Resuspend the cells in Isolation Buffer to 1×10^7 MNC per mL (or 1 mL for blood).
- For *positive isolation* or *depletion* of one cell type, add 25 µL pre-washed and resuspended magnetic beads. For *negative isolation* (removal of multiple cell types simultaneously) add 100 µL magnetic beads.
- Incubate for 20 min (positive isolation) or 30 min (depletion/negative isolation) at 2°C to 8°C with gentle tilting and rotation.
- Optional:* Add 1 mL Isolation Buffer to limit trapping of unbound cells.

8. Place the tube in a magnet for 2 min.
9. *Negative isolation/depletion:* Transfer the supernatant containing the unbound cells to a fresh tube for further experiments.
or
Positive isolation: While the tube is still in the magnet, carefully remove and discard the supernatant.
10. Remove the tube from the magnet and add 1 mL Isolation Buffer, pipet 2–3 times (or vortex 2–3 sec) and place the tube in a magnet for 2 min.
11. Repeat steps 9–10 at least twice to wash the cells and obtain high purity.
12. Resuspend the cell pellet in preferred cell medium.

Table 1: Volumes for indirect cell isolation.

Step	Step description	Volumes per 1 × 10 ⁷ MNC	Volumes per 2 × 10 ⁸ MNC
	Recommended tube size	5–7 mL	50 mL
	Recommended magnet	DynaMag™-5	DynaMag™-50
1	Primary biotinylated antibody	~10 µg	~200 µg
1	Cell volume (MNC/blood)	1 mL	20 mL
3*	Wash cells (Isolation Buffer)	~2 mL	~40 mL
4	Resuspend cells	1 mL	20 mL
5**(*)	Add Dynabeads® magnetic beads (positive isolation/depletion) Add Dynabeads® magnetic beads (negative isolation)	25 µL 100 µL	500 µL 2 mL
7*	Increase volume (Isolation Buffer)	~1 mL	~20 mL
10–12*	For positive isolation only: Wash cells (Buffer 1)	3 × ~1 mL	3 × 20 mL

* Adjust the Isolation Buffer volumes to fit to the tube you are using.

** If very high depletion efficiency is required or you are depleting many cells simultaneously, you might have to increase/ optimize the amount of Dynabeads® magnetic beads.

*** When incubating, tilt and rotate so the cells and beads are kept in the bottom of the tube. Do not perform end-over-end mixing if the volume is small relative to the tube size.

Isolate Cells – Direct Technique (Antibody-coating of Dynabeads® Magnetic Beads)

- See "General Guidelines" for recommendation of when to use the direct vs. indirect cell isolation technique.
- Use 0.5–1.5 µg of primary biotinylated antibody per 25 µL (1 × 10⁷) Dynabeads® magnetic beads. Titrate the primary antibody to optimize the amount used.

This protocol is based on 1 × 10⁷ MNC or 1 mL whole blood, but is directly scalable from 1 × 10⁷ to 4 × 10⁸ cells or 1–40 mL whole blood. When working with lower volumes than 1 × 10⁷ cells or 1 mL blood, use the same volumes as for 1 × 10⁷ cells or 1 mL blood. When working with larger volumes, scale up all reagent and volumes accordingly, as shown in Table 2.

1. Transfer 25 µL pre-washed and resuspended magnetic beads to a tube.
2. Add ~1 µg antibodies (titrate the antibody amount for your use).
3. Incubate for ≥30 minutes at room temperature with gentle tilting and rotation.
4. Place the tube in a magnet for 1 min and discard the supernatant.
5. Remove the tube from the magnet and add 2 mL Isolation Buffer.
6. Repeat steps 4–5 once to remove excess of antibodies.
7. Place the tube in the magnet for 1 min and discard the supernatant.
8. Add to the beads to 1 mL cell sample (blood or 1 × 10⁷ MNC) and resuspend.
9. Incubate for 20 min (positive isolation) or 30 min (depletion) at 2°C to 8°C with gentle tilting and rotation
10. *Optional:* Add 1 mL Isolation Buffer to limit trapping of unbound cells.
11. Place the tube in a magnet for 2 min.
12. *Depletion:* Transfer the supernatant containing the unbound cells to a fresh tube for further experiments.
or
Positive isolation: While the tube is still in the magnet, carefully remove and discard the supernatant.
13. Remove the tube from the magnet and add 1 mL Isolation Buffer, pipet 2–3 times (or vortex 2–3 sec) and place the tube in a magnet for 2 min.
14. Repeat steps 12–13 at least twice to wash the cells and to obtain high purity.
15. Resuspend the cell pellet in preferred cell medium.

Table 2: Volumes for direct cell isolation.

Step	Step description	Volumes per 1 × 10 ⁷ MNC	Volumes per 2 × 10 ⁸ MNC
	Recommended tube size	5–7 mL	50 mL
	Recommended magnet	DynaMag™-5	DynaMag™-50
1**	Dynabeads® magnetic beads	25 µL	500 µL
2**	Primary biotinylated antibody	~1 µg	~20 µg
4*	Isolation Buffer	1 mL	20 mL
5–6	Wash magnetic beads (Isolation Buffer)	2 x ~2 mL	2 x ~40 mL
8***	Cell volume	1 mL	20 mL
10*	<i>Optional:</i> Increase volume (Isolation Buffer)	~1 mL	~8 mL
13-14*	For positive isolation only: Wash the cells (Isolation Buffer)	3 × ~1 mL	3 × 20 mL

* Adjust the Buffer volumes to fit to the tube you are using.

** If the target cell population is high (e.g. >2.5 × 10⁶ target cells/mL), increase/double the amount of Dynabeads® magnetic beads and antibody.

*** When incubating, tilt and rotate so the cells and beads are kept in the bottom of the tube. Do not perform end-over-end mixing if the volume is small relative to the tube size.

Description of Materials

Dynabeads® Biotin Binder are uniform, superparamagnetic polystyrene beads (2.8 µm diameter) coated with with recombinant streptavidin. The streptavidin coated onto Dynabeads® magnetic beads will bind most biotinylated ligands. Unwanted binding of cells to streptavidin via lectin-like receptors or other adhesive receptors is avoided since the recombinant streptavidin contains neither sugar nor the RYD sequence.

Related Products

Product	Cat. no.
DynaMag™-5	12303D
DynaMag™-15	12301D
DynaMag™-50	12302D
Dynabeads® FlowComp™ Flexi	11061D
CELLlection™ Biotin Binder	11533D
HulaMixer® Sample Mixer	15920D

[REF] on labels is the symbol for catalog number.

Limited Use Label License

The purchase of this product conveys to the purchaser the limited, nontransferable right to use the purchased amount of the product only to perform internal research for the sole benefit of the purchaser. No right to resell this product or any of its components is conveyed expressly, by implication, or by estoppel. This product is for internal research purposes only and is not for use in commercial applications of any kind, including, without limitation, quality control and commercial services such as reporting the results of purchaser's activities for a fee or other form of consideration. For information on obtaining additional rights, please contact outlicensing@lifetech.com or Out Licensing, Life Technologies, 5791 Van Allen Way, Carlsbad, California 92008. Manufactured by Life Technologies AS, Norway. Life Technologies AS complies with the Quality System Standards ISO 9001:2008 and ISO 13485:2012.

Limited Product Warranty

Life Technologies Corporation and/or its affiliate(s) warrant their products as set forth in the Life Technologies' General Terms and Conditions of Sale found on Life Technologies' website at www.lifetechnologies.com/termsandconditions. If you have any questions, please contact Life Technologies at www.lifetechnologies.com/support.

SPEC-05724

©2013 Life Technologies Corporation. All rights reserved. The trademarks mentioned herein are the property of Life Technologies Corporation and/or its affiliate(s) or their respective owners.

DISCLAIMER: LIFE TECHNOLOGIES CORPORATION AND/OR ITS AFFILIATE(S) DISCLAIM ALL WARRANTIES WITH RESPECT TO THIS DOCUMENT, EXPRESSED OR IMPLIED, INCLUDING BUT NOT LIMITED TO THOSE OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE OR NON-INFRINGEMENT. TO THE EXTENT ALLOWED BY LAW, IN NO EVENT SHALL LIFE TECHNOLOGIES AND/OR ITS AFFILIATE(S) BE LIABLE, WHETHER IN CONTRACT, TORT, WARRANTY, OR UNDER ANY STATUTE OR ON ANY OTHER BASIS FOR SPECIAL, INCIDENTAL, INDIRECT, PUNITIVE, MULTIPLE OR CONSEQUENTIAL DAMAGES IN CONNECTION WITH OR ARISING FROM THIS DOCUMENT, INCLUDING BUT NOT LIMITED TO THE USE THEREOF.

For support visit www.lifetechnologies.com/support or email techsupport@lifetech.com

www.lifetechnologies.com

life
technologies™